

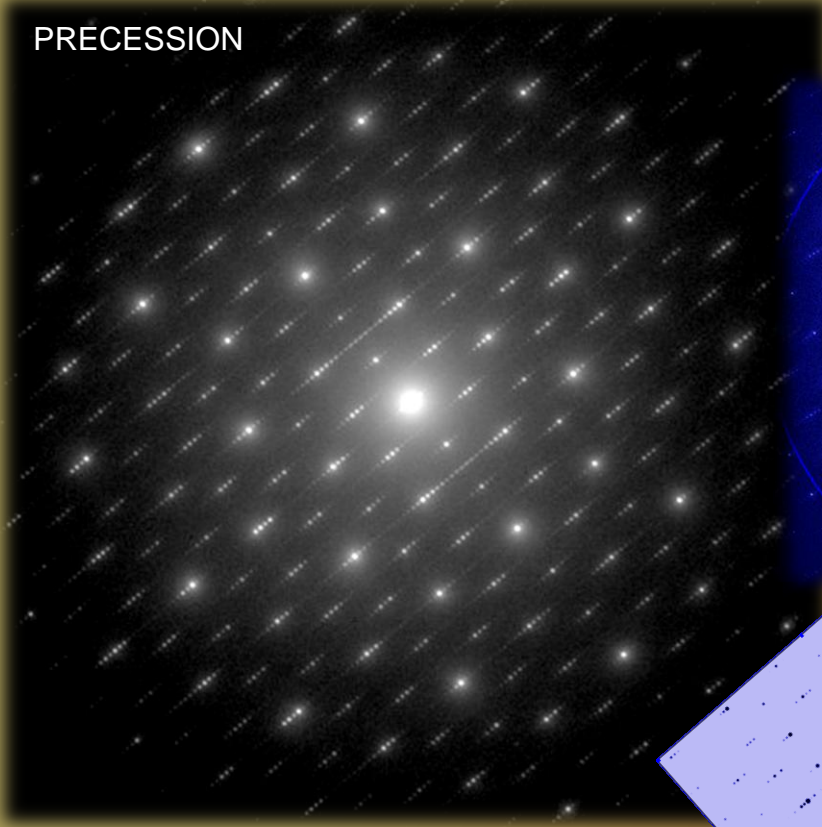



PRECESSION

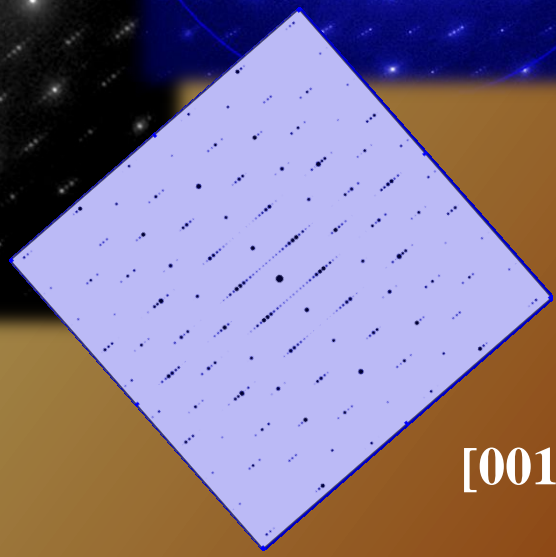

NO PRECESSION

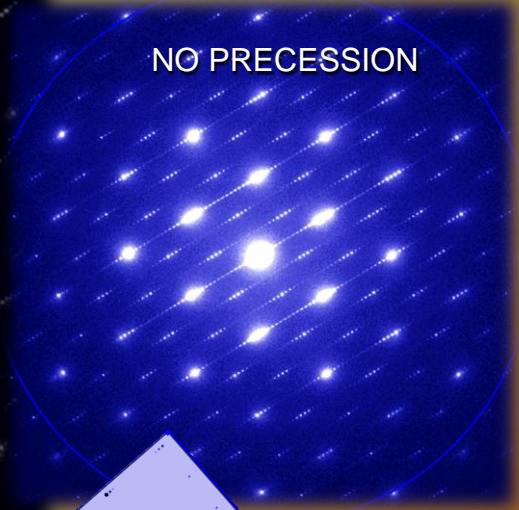

[001]

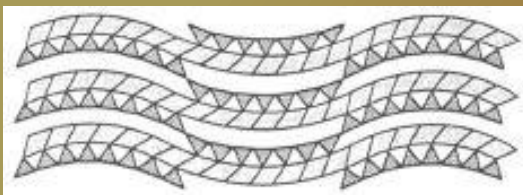

Antigorite is a serpentine mineral ($\sim\text{Mg}_3\text{Si}_2\text{O}_5(\text{OH})_4$) with a modulated, wavy structure. Using precession ED patterns are closer to kinematical (see simulation in the middle) allows to find the nature of the superstructure

“spinning star” precession device JEOL 2010 200kv, precession angle 1.5° , Gatan TV rate camera, FUJI image plates

COURTESY: Dr. G.Capitani Univ of Bari (Italia), Prof. M.Mellini Univ.of Siena (Italia)

